

Check veterinary drug & pesticides in honey with FartakTeb Rapid Tests Kit

- ✿ One Step Universal Preparation
- ✿ Rapid 10min Test
- ✿ No Instrument Required




Drug Residue Rapid Test Kit


Item #	Product - 48 tests/kit
100307	Tetracyclines Rapid Test Kit
100309	Quinolones Rapid Test Kit Kit
100314	Chloramphenicol Rapid Test Kit
100323	Tilmicosin & Tylosin, Macrolides Rapid Test Kit
100334	Carbendazim Rapid Test Kit
100335	Acetamiprid Rapid Test Kit
100336	Spectinomycin, Thiamphenicol & Quinolones, SQT Rapid Test Kit
100337	Tetracyclines, Macrolides & Erythromycin, TME Rapid Test Kit
100338	Carbendazim, Acetamiprid & Carbaryl, MDL Rapid Test Kit
100339	Imidacloprid, Cyromazine, Thiamethoxam, Carbofuran, DEMN Rapid Test Kit

We also have ELISA kit for honey testing. Contact for details.

Procedures for testing

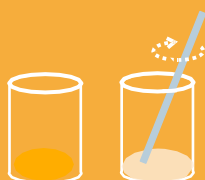

Take honey into a tube and mix with sample buffer

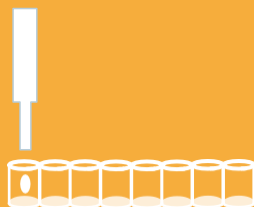

Add 200ul mixed sample into vial mix and incubate at 40°C
Wait for 5min



Insert test strip
Wait for 5min
at 40°C



Check the result

** this procedure applies for drug residue with MRL.

For not allowed drugs such as nitrofurantoin, please contact us



FARTAKTEB Co Ltd

Add: Vanak SQ., Gandhi Jonoubi ST., No. 146

Tel: +98-9122904915

Email: info@fartak-teb.com

Web: www.fartak-teb.ir