



Check the meat and shrimp with Antibiotic Residue Rapid Test Kit

- 🌸 One Step
- 🌸 Rapid 10min Test
- 🌸 No Instrument Required




Antibiotic Residue Rapid Test Kit

Item #	Product - 48 tests/kit	Sample
510001	Meat Beta Lactams Rapid Test Kit	Meat
510005	Meat Lincomycin Rapid Test Kit	Meat
510006	Meat Tetracyclines Rapid Test Kit	Meat
510008	Meat Sulfonamides Rapid Test Kit	Meat
510009	Meat Fluoroquinolones Rapid Test Kit	Meat
510014	Meat Chloramphenicol Rapid Test Kit	Meat
510023	Meat Macrolides Rapid Test Kit	Meat
510032	Meat Florfenicol Rapid Test Kit	Meat
510034	Meat Trimethoprim Rapid Test Kit	Meat
510035	Meat Olaquinox Rapid Test Kit	Meat
E4701	Nitrofurantoin AOZ/AMOZ/AHD/SEM 4-in-1 QuaTest	Shrimp
E1705	Chloramphenicol Rapid Test Kit	Shrimp

We also have ELISA Kits for meat and shrimp testing. Contact us for details.

Procedures for testing

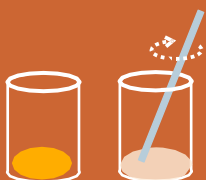

Homogenized meat into a tube and mix it with sample buffer



Centrifuge



Add 200ul transparent layer into vial mix and wait for 5min



Insert test strip
Wait for 5min
Check the result

**** this procedure applies for drug residue with MRL.**

For not allowed drugs such as CAP and nitrofurantoin, please contact us



FARTAKTEB Co Ltd

Add: Vanak SQ., Gandhi Jonoubi ST., No. 146

Tel: +98-9122904915

Email: info@fartak-teb.com

Web: www.fartak-teb.ir